

Innovative Products & Technologies
for the IVF/Fertility
Health Care Professionals



**THOMAS
MEDICAL**

a CRI® Company

Driven to innovate...focused on care™

6102 Victory Way Indianapolis, IN 46278 USA
800.556.0349 phone 888.411.3754 fax
sales@thomasmedical.com

www.thomasmedical.com

Tucker Embryo Transfer & Trial Catheter

Embryo Transfer, IVF

The Tucker Embryo Transfer Catheter is equipped with an easy loading ultra smooth tip

Specifications

FTR120 - Embryo Transfer Catheter, 20cm

FTR123 - Embryo Transfer Catheter, 23cm

FTR320 - Embryo Trial Transfer Catheter, 20cm

FTR323 - Embryo Trial Transfer Catheter, 23cm

Detail URM Tip




The ultrasonic refractory matte finish (URM) tip provides excellent unobstructed viewing of the embryo



20cm

23cm

URM Tip
Under Ultrasound

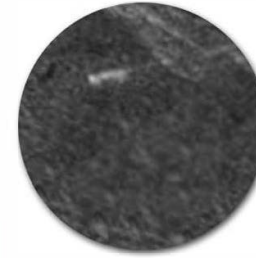

Superior Memory
Capabilities



Distinct Graduated
Markings



Unique Features of the Tucker Embryo Transfer Catheter

1. The unique design of the inner catheter prevents the Bernoulli effect negating turbulence and tumbling that can damage the embryo.

This design feature also negates the formation of bubbles in the column when loading and through the embryo transfer process

2. Crystal Clear inner catheter.
3. Catheter and (URM) echogenic tip enable clear view of the embryo at all times.
4. Outer sheath and inner catheter are designed with distinct graduated markings for ease of placement.
5. Outer sheath is equipped with superior memory characteristics affording excellent shapeability that is maintained through the procedure.
6. Inner catheter is designed with a universal luer fitting.
7. All catheters are one-cell mouse embryo bioassay (MEA) tested.




Fertility Essentials

Intrauterine Insemination, Diagnostic HSG Procedures, Endometrial Sampling

IUI Catheters

- TMI1192 - IUI Catheter, 5Fr, 15cm, Box of 25
- TMI1193 - Shapeable IUI Catheter, 5Fr, 15cm, Box of 25
- TMI1194 - IUI Catheter, 5Fr, 18cm, Box of 25
- TMI1195 - Shapeable IUI Catheter, 5Fr, 18cm, Box of 25
- TMI1295 - Clear Shapeable IUI Catheter, 5Fr, 18cm, Box of 25
- TMI1196 - IUI Catheter, 5Fr, 23cm, Box of 25
- TMI1197 - Shapeable IUI Catheter, 5Fr, 23cm, Box of 25
- TMI1297 - Clear Shapeable IUI Catheter, 5Fr, 23cm, Box of 25

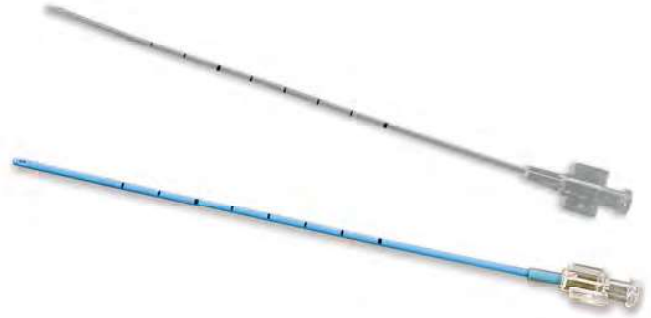


HS Catheters and HS Procedure Trays

- TMI1154 - Flexible HS Catheter, 5Fr, 28cm, Box of 10
- TMI1155 - Flexible HS Catheter, 7Fr, 28cm, Box of 10
- TMI1185 - Shapeable HS Catheter, 5Fr, 28cm, Box of 10
- TMI1187 - Shapeable HS Catheter, 7Fr, 28cm, Box of 10
- TMI1157 - HS Procedure Tray without Catheter, Box of 10
- TMI1160 - HS Procedure Tray with Flexible HS Catheter
5Fr, 28cm, Box of 10
- TMI1161 - HS Procedure Tray with Flexible HS Catheter
7Fr, 28cm, Box of 10
- TMI1198 - HS Procedure Tray with Shapeable HS Catheter
5Fr, 28cm, Box of 10
- TMI1199 - HS Procedure Tray with Shapeable HS Catheter
7Fr, 28cm, Box of 10



Endometrial Sampler

- TMI1175 Endometrial Sampler, Box of 25

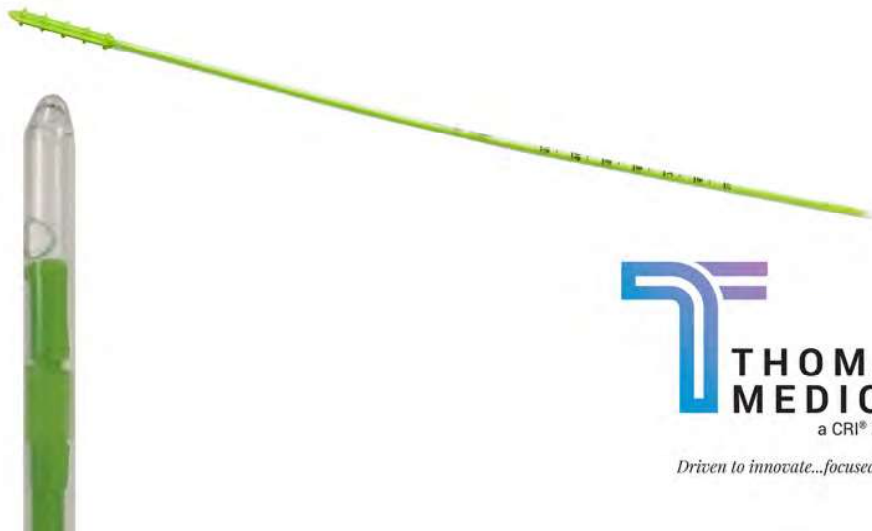

Driven to innovate...focused on care™

6102 Victory Way Indianapolis, IN 46278 USA
800.556.0349 phone 888.411.3754 fax
sales@thomasmedical.com

www.thomasmedical.com