


Designed to provide abundant trans-epithelial histology samples.



Soft ECC-S®

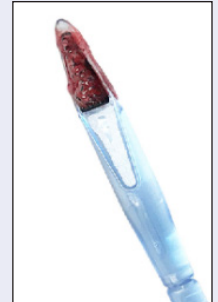

Soft ECC®

Examples of Soft ECC-S® and Soft ECC® Tissue Filled Tip.  
Goal: Entrap tissue into the hook array to ensure adequate samples.

## FEATURES

Ease of use: insert and rotate to gently scrape, not the old “insert and withdraw” sharp technique

Minimally invasive design

Abundant histological sample: simultaneous tissue collection and storage of the specimen for transport

Process as is customary for endocervical sample in the laboratory.

The Soft ECC-S® is a smaller pad size, designed for the shallow, short or stenotic cervix.

**WHAT IS KYLON®?** Kylon is a fabric with hooks that gently biopsies by compressing and frictionally lifting intact tissue fragments and simultaneously collects the specimen within the rows of hooks and fabric.

## TECHNICAL DATA

Length ..... 21cm

## ORDERING INFORMATION

HT-02-025 ..... Soft-ECC Endocervical Tissue Sampler

HT-03-025 ..... Soft-ECC-S Endocervical Tissue Sampler, Short

Manufactured by:



Compassionate Colposcopy®

TO ORDER:

www.thomasmedical.com | (800) 556-0349 | sales@thomasmedical.com  
6102 Victory Way Indianapolis, IN 46278

The Histologics Biopsy is a simple and gentle approach to perform a cervical or lower genital tract biopsy.



SoftBiopsy®

Example of SoftBiopsy® Tissue Filled Tip.  
Goal: Entrap tissue into the hook array to ensure adequate samples.

## FEATURES

Ease of use: Apply with moderate pressure to cervical or vaginal lesion and rotate to gently biopsy.

No pinching or cutting tissue

Minimally invasive design

Abundant histological sample: simultaneous tissue collection and storage of the specimen for transport

Process as is customary for exocervical biopsy sample in the laboratory.

**WHAT IS KYLON®?** Kylon is a fabric with hooks that gently biopsies by compressing and frictionally lifting intact tissue fragments and simultaneously collects the specimen within the rows of hooks and fabric.

## TECHNICAL DATA

Length ..... 20cm

## ORDERING INFORMATION

HT-01-025 ..... SoftBiopsy Exocervical Tissue Sampler

Manufactured by:



Compassionate Colposcopy®

TO ORDER:

www.thomasmedical.com | (800) 556-0349 | sales@thomasmedical.com  
6102 Victory Way Indianapolis, IN 46278