

Designed & manufactured to ensure a minimally invasive yet consistently effective procedure.

**FEATURES**

- Smooth tip reduces risk of trauma upon insertion and removal of catheter.
- Stable placement sheath decreases the need for a tenaculum.
- Soft, concentric balloon secures above the internal os preventing leakage of sterile saline or water-based contrast media.
- Shapeable HS catheter is equipped with an embedded stylet in cases when stenosis is present or when the vaginal tract and cervical canal are not aligned.
- Made with 20% Radiopaque material allowing catheter to be seen under xray.

Download IFU



**TECHNICAL DATA**

- French Size . . . . . 5Fr, 7Fr
- Length . . . . . 28cm
- Maximum Balloon Volume . . . . . 5Fr, 1.5ml/7Fr, 3ml
- Box Qty. . . . . 10 ea.
- Sterilization Method . . . . . R (Gamma)
- Standard Procedures. . . . . Hysterosalpingography, Saline Infusion Sonography, Sonohysterography

**ORDERING INFORMATION**

- TMI1154. . . . . 5Fr Flexible HS Catheter
- TMI1155. . . . . 7Fr Flexible HS Catheter
- TMI1185. . . . . 5Fr Shapeable HS Catheter
- TMI1187. . . . . 7Fr Shapeable HS Catheter

All Thomas Medical HS Catheters are PVC Free and Non-DEHP. Not made with natural rubber latex.

# SPECIALTY HS CATHETERS




## STEEN OPEN-TIP HS CATHETER

The Steen Open-Tip HS catheter utilizes a balloon closer to the distal tip in order to provide a clearer view of the uterus while performing a hysterosalpingography or sonohysterography procedure. The open distal end allows a forward injection into the uterus.

### ORDERING INFORMATION

- TMI1188..... 5Fr Steen Open-Tip HS Catheter
- TMI1189..... 7Fr Steen Open-Tip HS Catheter

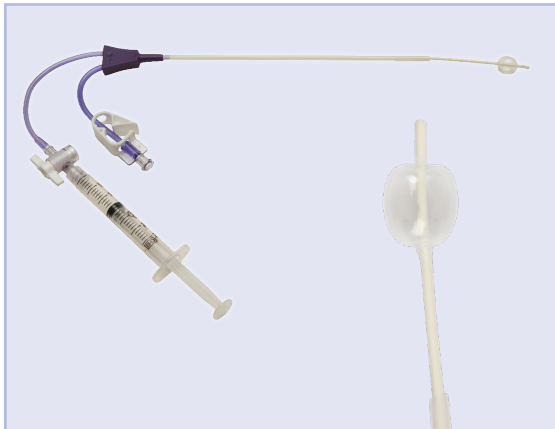


## MILLER ADVANCE HS CATHETER

The Miller Advance Catheter was designed with features for improved patient comfort with a soft catheter tip and slender catheter shaft that reduces the need for dilation.

### ORDERING INFORMATION

- TMI2000..... 5Fr, 21cm Miller Advance HS Catheter
- TMI2001..... 5Fr, 25cm Miller Advance HS Catheter
- TMI2002..... 5Fr, 28cm Miller Advance HS Catheter

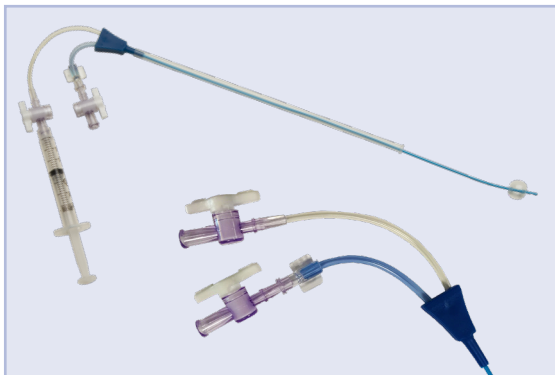

## TWO-TAP HS CATHETER

The Two-Tap HS catheter utilizes a second stopcock instead of a pinch clamp for easier manipulation when releasing the saline or water-based contrast media from the uterus.

### ORDERING INFORMATION

- TMI1181..... 5Fr Two-Tap HS Catheter
- TMI1182..... 7Fr Two-Tap HS Catheter
- TMI1183..... 5Fr Shapeable Two-Tap HS Catheter
- TMI1184..... 7Fr Shapeable Two-Tap HS Catheter

### TECHNICAL DATA

| CATHETER TYPE              | FRENCH SIZE(S) | LENGTH (CM) | MAXIMUM BALLOON VOLUME (ml) | BOX QTY | STERILIZATION METHOD |
|----------------------------|----------------|-------------|-----------------------------|---------|----------------------|
| FLEXIBLE HS CATHETER       | 5FR/7FR        | 28          | 5Fr, 1.5ml/7Fr, 3ml         | 10      | R (GAMMA)            |
| SHAPEABLE HS CATHETER      | 5FR/7FR        | 28          | 5Fr, 1.5ml/7Fr, 3ml         | 10      | R (GAMMA)            |
| STEEN OPEN-TIP HS CATHETER | 5FR/7FR        | 28          | 5Fr, 1.5ml/7Fr, 3ml         | 10      | R (GAMMA)            |
| MILLER ADVANCE HS CATHETER | 5FR            | 21, 25, 28  | 5Fr, 1.5ml                  | 10      | R (GAMMA)            |
| TWO-TAP HS CATHETER        | 5FR/7FR        | 28          | 5Fr, 1.5ml/7Fr, 3ml         | 10      | R (GAMMA)            |



Thomas Medical is an industry leading manufacturer of single-use women's health products.

www.thomasmedical.com | (800) 556-0349 | sales@thomasmedical.com  
6102 Victory Way Indianapolis, IN 46278 USA

